


Observer 9

Patient Monitor



www.dameca.com

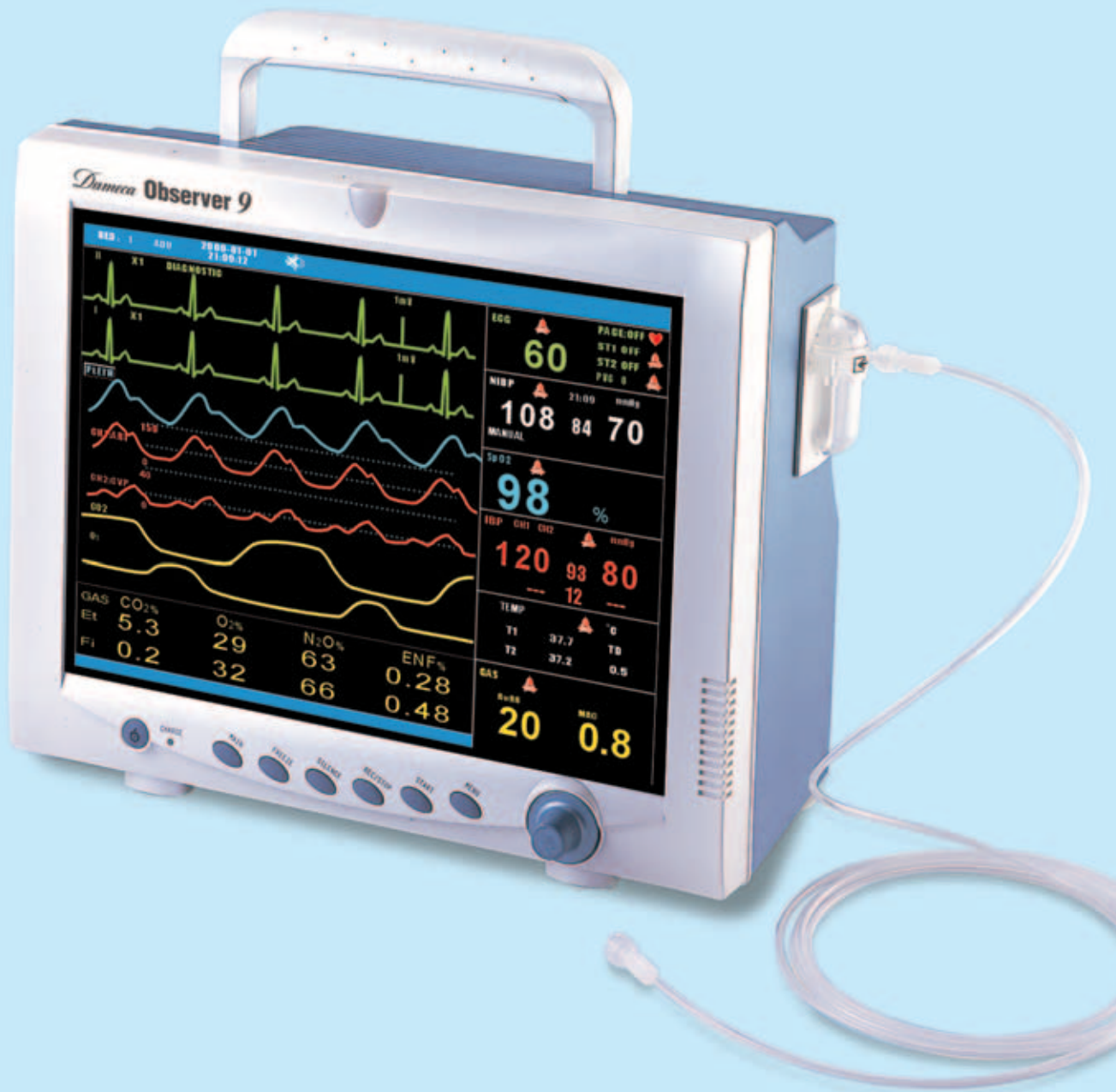

Dameca



Patient Monitor

The Observer patient monitors provide all essential anaesthesia monitoring features you need, whenever perioperative care is provided. Compact, easy to use, and the perfect partner for your anaesthesia workstation. Suitable for the full range of patients from neonatal to adult.

- Large 12.1" colour TFT display with maximum 8 waveforms
- Arrhythmia analysis, multi lead ST segment analysis and pacemaker detection
- Oxy CRG software available, optional
- 40 seconds of full-disclosure waveform review
- Drug dose calculation
- Networking solution with bed to bed view
- Multi language available *

Parameters: ECG, RESP, NIBP, SpO₂, 2 x TEMP, 2 x IBP, CO, CO₂ and Multi-Gas/O₂.



Multi-Gas/O₂ Measuring

Adopting infrared absorption technique to measure the concentration of Et and Fi values of CO₂, N₂O and five other anaesthetic agents including Des, Sev, Hal, Enf and Iso. O₂ measurement also available.



Large Colour Display

12.1" colour TFT display with 8 waveforms and parameter values in large size, featured by high resolution and particular design.



OxyCRG

Display of the interactive relation between heart rate, respiration and Oxygen parameter, judging the respiration and circulation function of neonate patients.

Observer 9



Trend Graph

Comprehensive review of monitoring information obtaining access to 72-hour trend graph of all parameters and 40 seconds full-disclosure waveforms.



Arrhythmia Review

Display of real-time waveforms and vital data produced prior to, during and after the occurrence of Arrhythmia events.



Rotary Knob

Ensuring accurate and prompt selection.



Multi-Gas/O₂ module

Fast response time. Automatic identification of anaesthetic agents.



Dameca

	3-5-lead ECG	ST Analysis	Imp. Resp.	SpO ₂ /Pleth	Temp x 2	NIBP	Oxy CRG	Network	CO ₂	O ₂	N ₂ O	Agents	Agent ID	IBP x 2/C.O.	Masimo SpO ₂	Recorder	Start-up Kit
Observer 9-1	●	●	●	●	●	●	●	●									●
Observer 9-2	●	●	●	●	●	●	●	●	●								●
Observer 9-3	●	●	●	●	●	●	●	●	●	●	●	●	●				●
Option A														●			
Option B															●		
Option C																●	

* **Languages Available:** Danish, English, French, German, Italian, Spanish, Portuguese and Polish.

Technical Data: Please see enclosed data sheet.



Dameca a/s - established in 1947- manufacturer of anaesthesia machines and supplier of patient monitors for the Danish and world markets. Dameca is centrally located 12 km from Copenhagen centre and 15 minutes drive from Copenhagen Airport.

We are constantly developing our products and consequently reserve the right to make alterations in design and equipment without prior notice. 04/05

Dameca a/s · Islevdalvej 211 · DK-2610 Roedovre · Copenhagen · Denmark

Phone: +45 4450 9990 · Fax: +45 4450 9999 · E-mail: info@dameca.com · www.dameca.com